


ADB-4en-PINACA

Sample Type: **Seized Material**

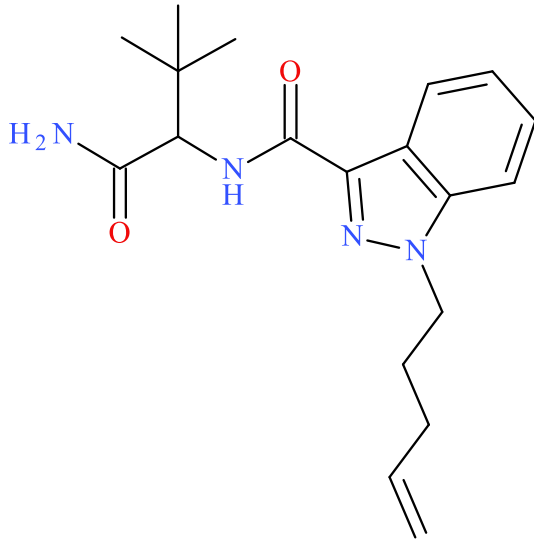

Latest Revision: **March 4, 2021**

Date Received: **January 8, 2021**

Date of Report: **March 4, 2021**

1. GENERAL INFORMATION

IUPAC Name:	N-(1-carbamoyl-2,2-dimethyl-propyl)-1-pent-4-enyl-indazole-3-carboxamide
InChI String:	InChI=1S/C19H26N4O2/c1-5-6-9-12-23-14-11-8-7-10-13(14)15(22-23)18(25)21-16(17(20)24)19(2,3)4/h5,7-8,10-11,16H,1,6,9,12H2,2-4H3,(H2,20,24)(H,21,25)
CFR:	Not Scheduled (03/2021)
CAS#	Not Available
Synonyms:	ADMB-4en-PINACA, ADB-PENINACA
Source:	NMS Labs – Criminalistic Laboratory
Appearance:	Plant-Like Material

Important Note: All identifications were made based on evaluation of analytical data (GC-MS and LC-QTOF-MS) in comparison to analysis of acquired reference material.

Prepared By: Alex J. Krotulski, PhD; Kristin Kendrick, MSFS; Sarah A. Shuda, MFSF, F-ABC; Melissa F. Fogarty, MSFS, D-ABFT-FT; Sarah E. Decker, BA; and Barry K. Logan, PhD, F-ABFT

2. CHEMICAL AND PHYSICAL DATA

2.1 CHEMICAL DATA

Form	Chemical Formula	Molecular Weight	Molecular Ion [M ⁺]	Exact Mass [M+H] ⁺
Base	C ₁₉ H ₂₆ N ₄ O ₂	342.4	342	343.2129

3. BRIEF DESCRIPTION

ADB-4en-PINACA is classified as a synthetic cannabinoid. Synthetic cannabinoids have been reported to cause psychoactive effects similar to delta-9-tetrahydrocannabinol (THC). Synthetic cannabinoids have caused adverse events, including deaths, as described in the literature. MDMB-4en-PINACA and ADB-PINACA are structurally similar synthetic cannabinoids. ADB-PINACA is a Schedule I substance in the United States; ADB-4en-PINACA and MDMB-4en-PINACA are not explicitly scheduled.

4. ADDITIONAL RESOURCES

<https://www.caymanchem.com/product/33205/adb-4en-pinaca>

5. QUALITATIVE DATA

5.1 GAS CHROMATOGRAPHY MASS SPECTROMETRY (GC-MS)

Testing Performed At: NMS Labs (Willow Grove, PA)

Sample Preparation: Acid/Base extraction

Instrument: Agilent 5975 Series GC/MSD System

Column: Zebron™ Inferno™ ZB-35HT (15 m x 250 μm x 0.25 μm)

Carrier Gas: Helium (Flow: 1 mL/min)

Temperatures: Injection Port: 265 °C
Transfer Line: 300 °C
MS Source: 230 °C
MS Quad: 150 °C

Oven Program: 60 °C for 0.5 min, 35 °C/min to 340 °C for 6.5 min

Injection Parameters: Injection Type: Splitless

Injection Volume: 1 µL

MS Parameters: Mass Scan Range: 40-550 m/z


Threshold: 250

Retention Time: 7.99 min

Standard Comparison: Reference material for ADB-4en-PINACA (Batch: 0606857-1) was purchased from Cayman Chemical (Ann Arbor, MI, USA). Analysis of this standard resulted in positive identification of the analyte in the exhibit as ADB-4en-PINACA based on retention time (7.96 min) and mass spectral data.



<https://www.caymanchem.com/product/33205/adb-4en-pinaca>

Chromatogram: ADB-4en-PINACA

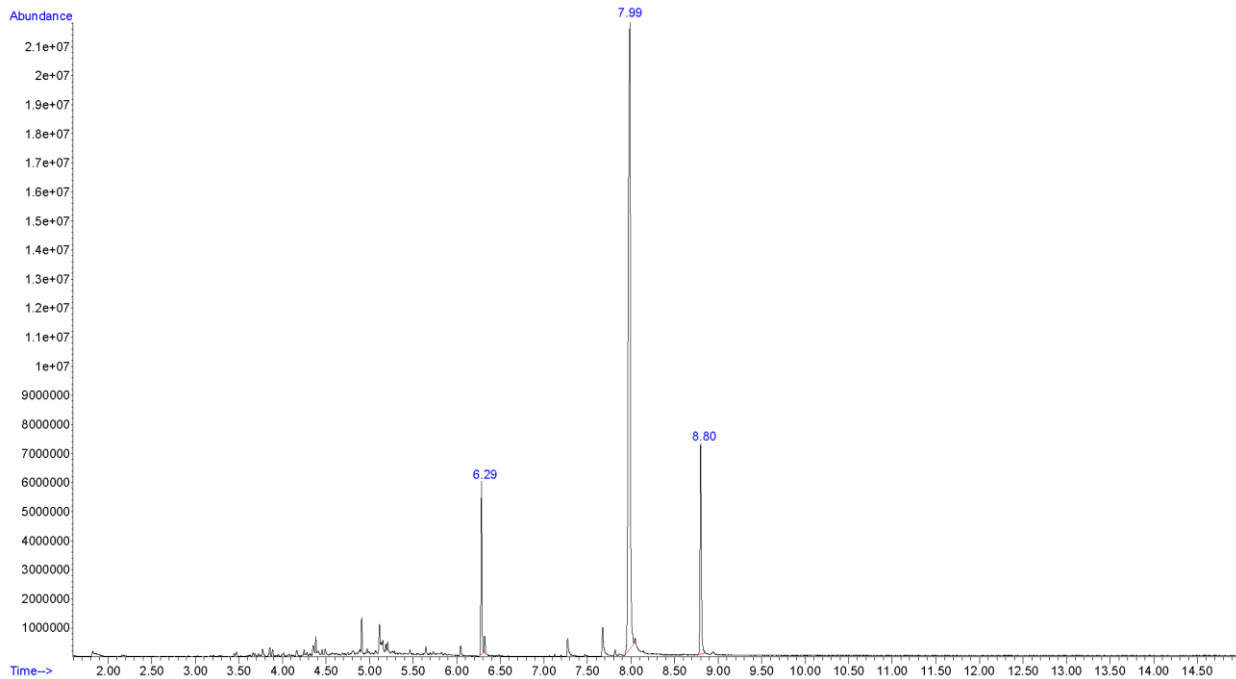

Additional peaks present in chromatogram: internal standard (6.29 min) and 4F-ABINACA (8.80 min)


EI (70 eV) Mass Spectrum (Top) and 10x (Bottom): ADB-4en-PINACA

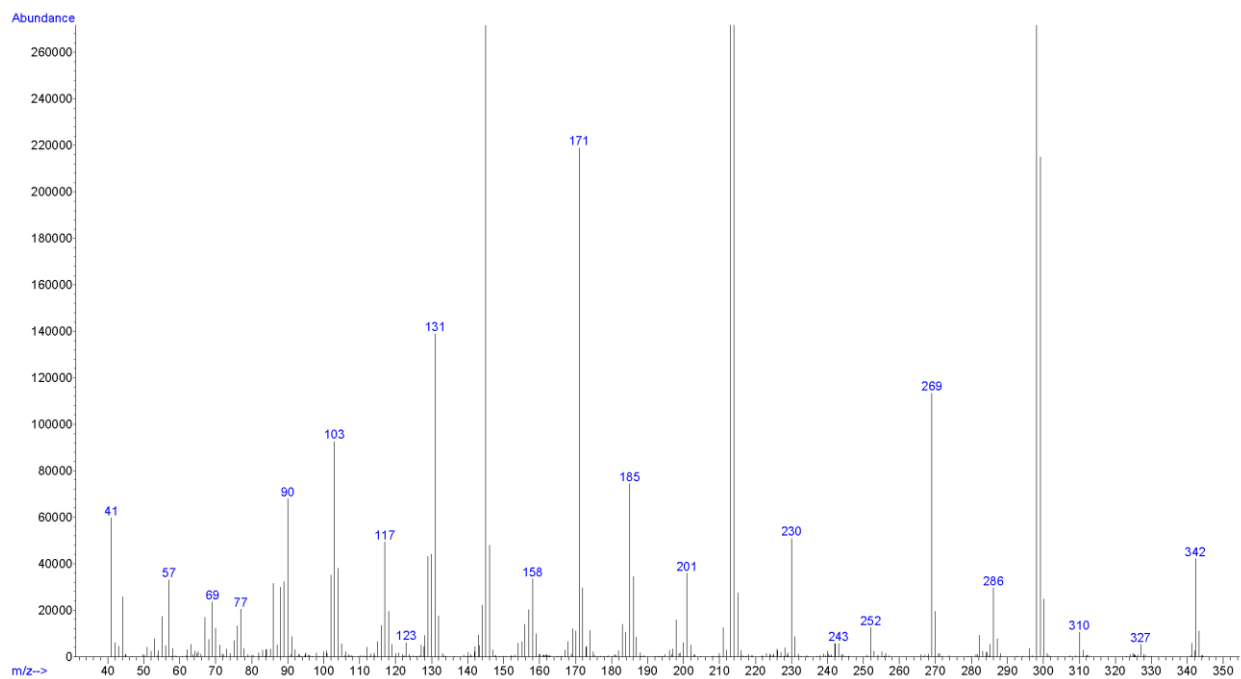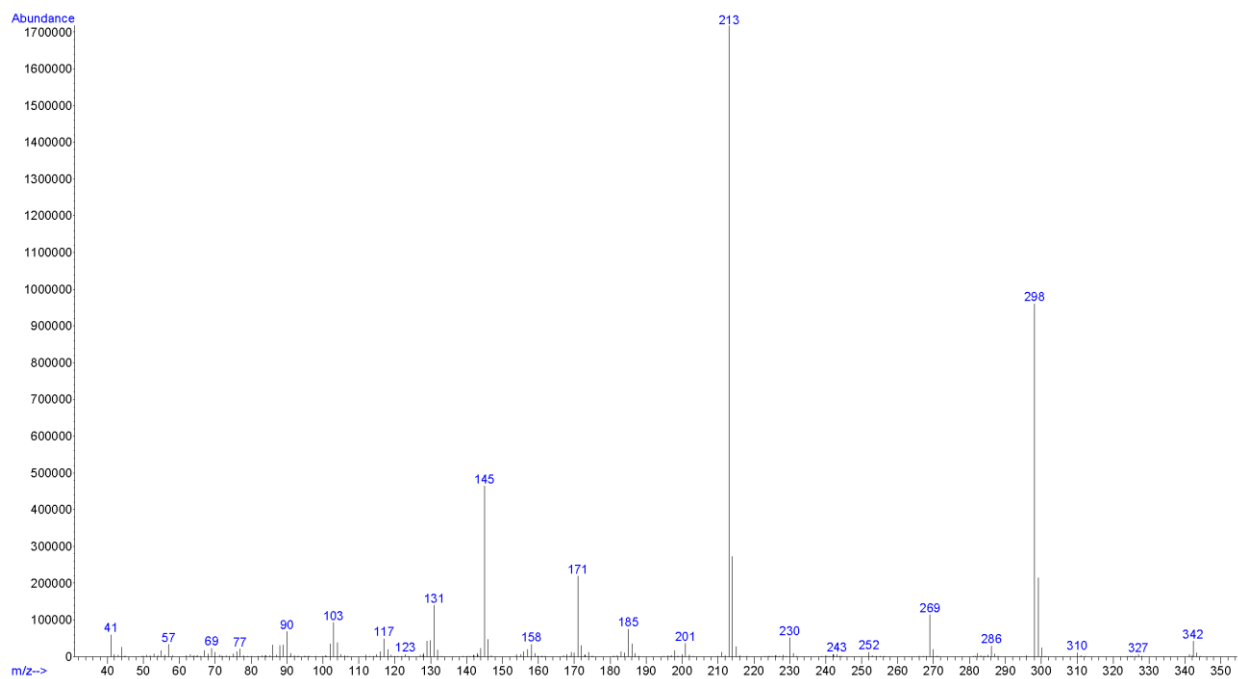

5.2 LIQUID CHROMATOGRAPHY QUADRUPOLE TIME OF FLIGHT MASS SPECTROMETRY (LC-QTOF)



Testing Performed At:	The Center for Forensic Science Research and Education at the Fredric Rieders Family Foundation (Willow Grove, PA)
Sample Preparation:	1:100 dilution of acid/base extract in mobile phase
Instrument:	Sciex TripleTOF® 5600+, Shimadzu Nexera XR UHPLC
Column:	Phenomenex® Kinetex C18 (50 mm x 3.0 mm, 2.6 µm)
Mobile Phase:	A: Ammonium formate (10 mM, pH 3.0) B: Methanol/acetonitrile (50:50) Flow rate: 0.4 mL/min
Gradient:	Initial: 95A:5B; 5A:95B over 13 min; 95A:5B at 15.5 min
Temperatures:	Autosampler: 15 °C Column Oven: 30 °C Source Heater: 600 °C
Injection Parameters:	Injection Volume: 10 µL
QTOF Parameters:	TOF MS Scan Range: 100-510 Da Precursor Isolation: SWATH® acquisition (27 windows) Fragmentation: Collision Energy Spread (35±15 eV) MS/MS Scan Range: 50-510 Da
Retention Time:	8.99 min
Standard Comparison:	Reference material for ADB-4en-PINACA (Batch: 0606857-1) was purchased from Cayman Chemical (Ann Arbor, MI, USA). Analysis of this standard resulted in positive identification of the analyte in the exhibit as ADB-4en-PINACA based on retention time (8.99 min) and mass spectral data. (https://www.caymanchem.com/product/33205/adb-4en-pinaca)

Chromatogram: ADB-4en-PINACA

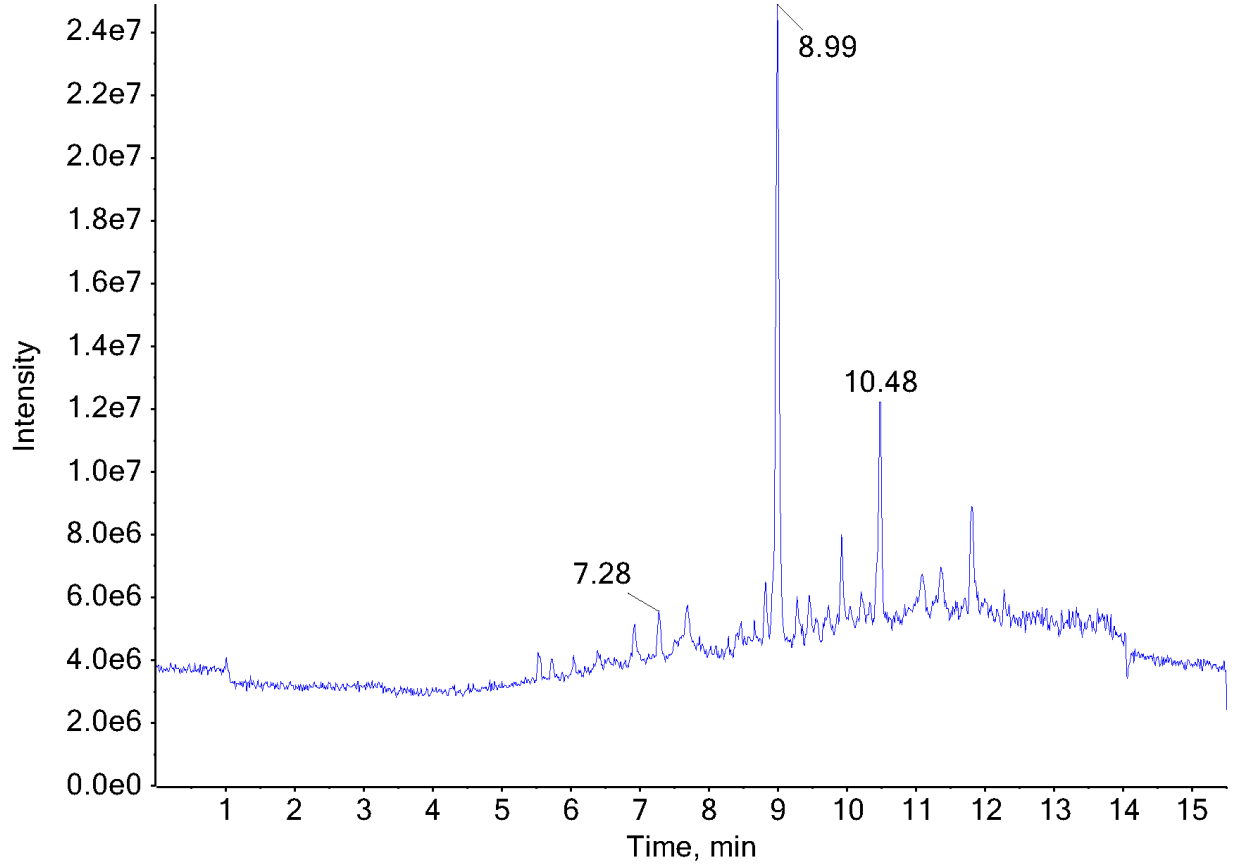

Additional peaks present in chromatogram: internal standard (7.28 min) and 4F-ABINACA (10.48 min)

TOF MS (Top) and MS/MS (Bottom) Spectra: ADB-4en-PINACA

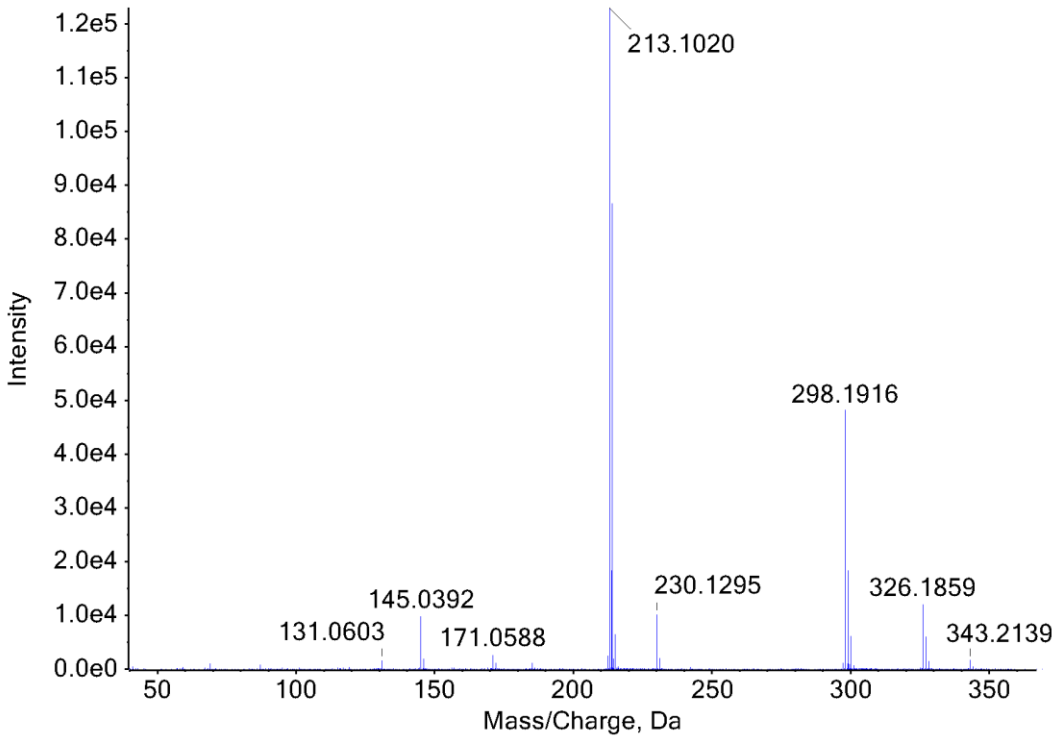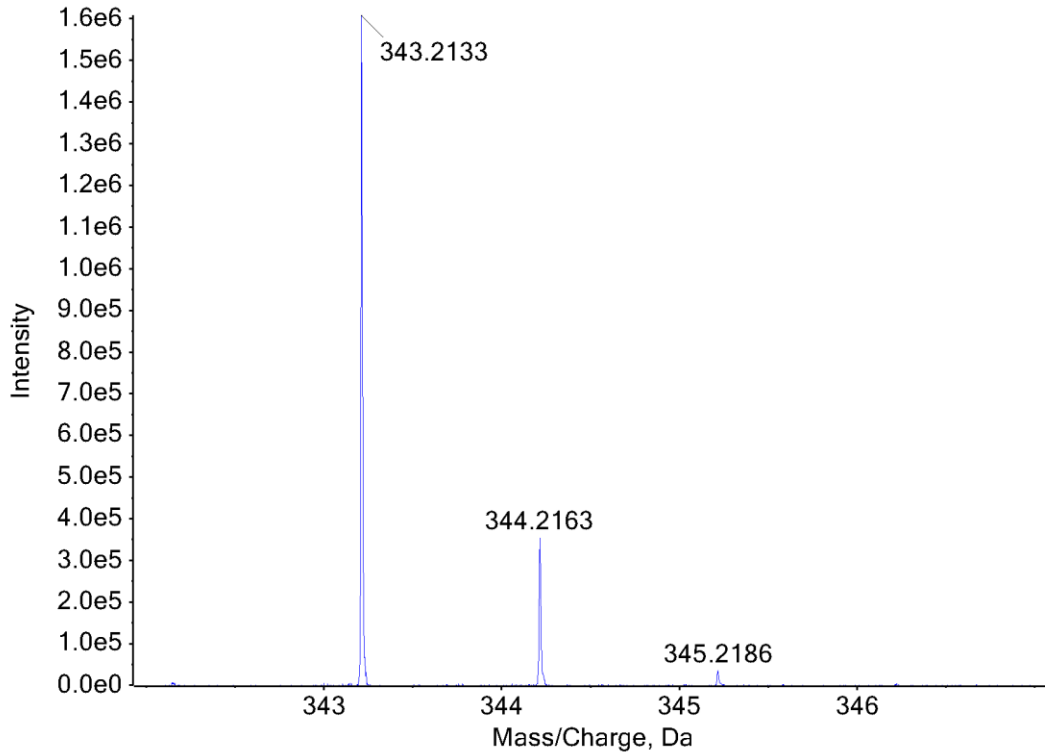

6. FUNDING

NPS Discovery at the Center for Forensic Science Research and Education (CFSRE) is supported in part by the National Institute of Justice, Office of Justice Programs, U.S. Department of Justice (Award Number 2020-DQ-BX-0007, “Real-Time Sample-Mining and Data-Mining Approaches for the Discovery of Novel Psychoactive Substances (NPS)”). The opinions, findings, conclusions and/or recommendations expressed in this publication are those of the author(s) and do not necessarily reflect those of the Department of Justice.